

# **CUBE 30 ATV**

A NEW APPROACH TO NON-INVASIVE VENTILATION

## **TECHNICAL SPECIFICATIONS**

# Ventilation modes CPAP

PSV S / BiLevel S PCV / BiLevel T PSV / BiLevel ST APCV

Pressure range

IPAP 5-30 / EPAP 3-25 cmH20 CPAP 3-20 cmH20

Pressure support 0-25 cmH20

Target Volume: OFF, 200 - 2.000ml

**Target Minute Volume:** OFF, 2.000 - 15.000ml

**Respiratory rate** 4-30 bpm

Inspiration time 0.3 - 5.0 seconds

**Inspiration Trigger**Five Level settings

Five Level settings

**Expiration Trigger** Auto, Manual settings 70%-30%

**Trigger Lock** 0.3 - 8.0 seconds

# Alarms

11 System/Therapy alarms

Power supply Internal Switch Mode Power Supply 100-240 V, 50/60 Hz

**Dimensions**Size (W x H x D)
168 x 108 x 182mm

Weight 1,6 Kg

Made in Germany

# Maximum Performance, minimum footprint & Made in Germany!

The NEW Cube 30 ATV from Drive DeVilbiss International is an innovative non-invasive volume supporting ventilator. The Cube 30 ATV incorporates intuitive state of the art functionality allowing safe and easy application of the device in a homecare or clinical environment. With its unique trigger technology and target minute volume algorithm, the CUBE 30 ATV will optimise therapy to meet individual patient requirements.

#### **DESIGN**

■ The device includes a large 4.3" coloured screen display incorporating easy to use tile design menu screens for clear convenient user navigation.






# TECHNOLOGY

- High sensitive flow based trigger detection with trigger failure auto-detect function.
- Contemporary volume assured pressure support which can be set based on breath-by-breath target volume or target minute volume.
- Extensive device alarm management, including on screen display alarm protocol.
- Swiss made turbine blower promising high quality performance for extended, reliable operation.

## **INFORMATIVE**

- Showing clear simple text information the device displays corresponding data settings including user notifications for incorrect parameter settings.
- Three different menu screen colours identifying Patient and Clinical menu modes.



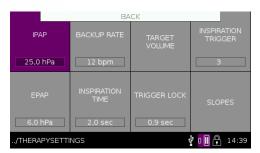

#### **VARIOUS ADJUSTMENT OPTIONS**

- Two independent adjustable ventilation settings allowing patients to change therapy mode for daytime and night time use.
- The Mask Check function ensures a correct fitting mask even before therapy starts.
  On screen values identify leakages in order for the mask to be adjusted and corrected.





# **OPTIONAL HUMIDIFIER**

User friendly and high humidity performance for maximum comfort during therapy.





The integrated humidifier connects underneath the Cube 30 ATV without increasing the foot print of the device and therefore saving valuable space.

**Drive DeVilbiss International** 



